

**DRG<sup>®</sup> Insulin (Porcine) ELISA (EIA-4747)****As of 10 Sept. 2007****For Veterinary Use Only****INTENDED USE**

The Insulin (Porcine) ELISA provides a method for the quantitative determination of insulin in porcine serum and plasma.

**SUMMARY AND EXPLANATION OF THE TEST**

Insulin is the principal hormone responsible for the control of glucose metabolism. It is synthesised in the  $\beta$ -cells of the islets of Langerhans as the precursor, proinsulin, which is processed to form C-peptide and insulin. Both are secreted in equimolar amounts into the portal circulation. The mature insulin molecule comprises two polypeptide chains, the A chain and the B chain. The two chains are linked together by two inter-chain disulphide bridges. There is also an intra-chain disulphide bridge in the A chain.

Secretion of insulin is mainly controlled by plasma glucose concentration, and the hormone has a number of important metabolic actions. Its principal function is to control the uptake and utilisation of glucose in peripheral tissues via the glucose transporter. This and other hypoglycaemic activities, such as the inhibition of hepatic gluconeogenesis and glycogenolysis are counteracted by the hyperglycaemic hormones including glucagon, epinephrine (adrenaline), growth hormone and cortisol.

Islet cell transplantation has long been considered as a potential cure for type 1 diabetes. The shortage of human donors and difficulty in isolating purified islets from adult human pancreas has drawn the attention to porcine pancreatic islets (1-4). Pig and human insulin are structurally similar, and the regulation of insulin secretion in pigs resembles humans. Parameters such as insulin secretion are used to investigate the viability of the transplanted islets (4-6).

**PRINCIPLE OF THE PROCEDURE**

The Insulin (Porcine) ELISA is a solid phase two-site enzyme immunoassay. It is based on the direct sandwich technique in which two monoclonal antibodies are directed against separate antigenic determinants on the insulin molecule.

During incubation, insulin in the sample reacts with peroxidase-conjugated anti-insulin antibodies and anti-insulin antibodies bound to the microtitration well. After a simple washing step that removes unbound enzyme labelled antibody, the bound conjugate is detected by reaction with 3,3'-5,5'-tetramethylbenzidine (TMB). The reaction is stopped by the addition of acid, giving a colorimetric endpoint that can be read spectrophotometrically.

**WARNINGS AND PRECAUTIONS**

- The content of this kit and their residues must not be allowed to come into contact with ruminating animal or swine.
- The Stop Solution in this kit contains 0.5 M H<sub>2</sub>SO<sub>4</sub>. Follow routine precautions for handling hazardous chemicals.
- All patient samples should be handled as capable of transmitting infections.

**MATERIAL REQUIRED BUT NOT PROVIDED**

- Pipettes for 25, 50, 100, 200 and 1000  $\mu$ l (repeat pipettes preferred for addition of enzyme conjugate solution, TMB Substrate and Stop Solution)
- Beakers and cylinders for reagent preparation

---

**DRG International Inc., USA****Fax: (908) 233-0758 • E-mail: [corp@drg-international.com](mailto:corp@drg-international.com) • Web: [www.drg-international.com](http://www.drg-international.com)**

## DRG<sup>®</sup> Insulin (Porcine) ELISA (EIA-4747)

As of 10 Sept. 2007

For Veterinary Use Only

- Redistilled water
- Microplate reader (450 nm filter)
- Plate shaker (The recommended velocity is 700-900 cycles per minute, orbital movement)
- Microplate washing device

### REAGENTS

Each Insulin (Porcine) ELISA kit contains reagents for 96 wells, sufficient for 42 samples and one standard curve in duplicate. For larger series of assays, use pooled reagents from packages bearing identical lot numbers.

The expiry date for the complete kit is stated on the outer label.

The recommended storage temperature is 2-8°C.

<b>Coated Plate</b> (Mouse monoclonal anti-insulin)	1 plate, 8-well strips	96 wells	Ready for use For unused microtitration strips, reseal the bag using adhesive tape and store at 2-8°C for 2 months.
<b>Standard 1, 2, 3, 4, 5</b> (porcine insulin) Concentration stated on vial label.	5 vials	1000 µl	Ready for use
<b>Standard 0</b> Color coded yellow	1 vial	5 ml	Ready for use
<b>Enzyme Conjugate 11X</b> (Peroxidase conjugated mouse monoclonal anti-insulin)	1 vial	1.3 ml	<i>Preparation, see below</i>
<b>Enzyme Conjugate Buffer</b> Color coded blue	1 vial	13 ml	Ready for use
<b>Wash Buffer 21X</b> Storage after dilution: 2-8°C for 2 months	1 bottle	40 ml	<i>Dilute</i> with 800 ml redistilled water to make wash buffer.
<b>TMB Substrate Solution</b> (TMB) Colorless solution <i>Note! Light sensitive!</i>	1 bottle	22 ml	Ready for use
<b>Stop Solution</b> 0.5 M H <sub>2</sub> SO <sub>4</sub>	1 vial	7 ml	Ready for use

### 1.1 Preparation of enzyme conjugate solution

Prepare the needed volume of enzyme conjugate solution by dilution of Enzyme Conjugate 11X, (1+10) in Enzyme Conjugate Buffer according to the table below. Mix gently.

**DRG<sup>®</sup> Insulin (Porcine) ELISA (EIA-4747)****As of 10 Sept. 2007****For Veterinary Use Only**

When preparing enzyme conjugate solution for the whole plate or if the reagents are to be used within 2 weeks, pour all of the Enzyme Conjugate Buffer into the Enzyme Conjugate 11X vial.

Number of strips	Enzyme Conjugate 11X	Enzyme Conjugate Buffer
12 strips	1 vial	1 vial
8 strips	700 µl	7 ml
6 strips	500 µl	5 ml
4 strips	400 µl	4 ml

Storage after dilution: 2-8°C for 2 weeks.

**SPECIMEN COLLECTION AND HANDLING****Serum**

Collect blood by venipuncture, allow to clot and separate the serum by centrifugation.

Samples can be stored at 2-8°C up to 24 hours. For longer periods, store samples at -20°C. Avoid repeated freezing and thawing.

**Plasma**

Collect blood by venipuncture into tubes containing heparin, citrate or EDTA as anticoagulant, and separate the plasma fraction.

Samples can be stored at 2-8°C up to 24 hours. For longer periods store samples at -20°C. Avoid repeated freezing and thawing.

**1.2 PREPARATION OF SAMPLES**

No dilution is normally required for serum or plasma.

All samples containing porcine insulin above the highest standard should be diluted with Standard 0 or with Diabetes Sample Buffer [REF] 10-1195-01.

**TEST PROCEDURE**

All reagents and samples must be brought to room temperature before use. Perform each determination in duplicate for standards and unknowns. Prepare a standard curve for each assay run.

1. Prepare enzyme conjugate solution (according to the table on previous page) and wash buffer.
2. Prepare sufficient microplate wells to accommodate Standards and samples in duplicate.
3. Pipette 25 µl each of Standards and samples into appropriate wells.
4. Add 100 µl of enzyme conjugate solution into each well.
5. Incubate on a plate shaker (700-900 rpm) for 2 hours at room temperature (18-25°C).

## DRG<sup>®</sup> Insulin (Porcine) ELISA (EIA-4747)

As of 10 Sept. 2007

For Veterinary Use Only

6. Wash plate 6 times with 700 µl wash buffer per well with an automatic washer. \*  
After final wash, invert and tap the plate against absorbent paper.
7. Add 200 µl TMB Substrate Solution into each well.
8. Incubate for 15 minutes at room temperature (18-25°C).
9. Add 50 µl Stop Solution to each well.  
Place the plate on the shaker for approximately 5 seconds to ensure mixing.
10. Read optical density at 450 nm and calculate results.  
Read within 30 minutes.

\* The plate can be washed with an automatic washer or manually. When washing with an automatic washer please use the overflow function. If there is no overflow function available on the automatic washer please wash manually.

Manual wash can be done either with a pipette or a squirt bottle: Aspirate the reaction volume and add 400 µl wash buffer to each well with a pipette or fill the wells completely by spraying wash buffer into the wells with a squirt bottle. Aspirate completely and repeat 5 times. The overflow is not a problem rather an advantage.

Note! To prevent contamination between the conjugate and substrate, separate pipettes are recommended.

### INTERNAL QUALITY CONTROL

Commercial control and/or internal serum pools with low, intermediate and high porcine insulin concentrations should routinely be assayed as unknowns, and results charted from day to day.

It is good laboratory practice to record the following data for each assay: kit lot number; reconstitution dates of kit components; OD values for the blank, Standards and controls.

### CALCULATION OF RESULTS

#### Computerized calculation

The concentration of porcine insulin is obtained by computerized data reduction of the absorbance for the Standards, except Standard 0, versus the concentration using cubic spline regression.

#### Manual calculation

1. Plot the absorbance values obtained for the Standards, except Standard 0, against the porcine insulin concentration on a lin-lin paper and construct a standard curve.
2. Read the concentration of the unknown samples from the standard curve.

## DRG<sup>®</sup> Insulin (Porcine) ELISA (EIA-4747)

As of 10 Sept. 2007

For Veterinary Use Only

### Example of results

Wells	Identity	A <sub>450</sub>	Mean conc. µg/l
1A-B	Standard 0	0.054/0.053	
1C-D	Standard 1 (0.02 µg/l)*	0.090/0.086	
1E-F	Standard 2 (0.05 µg/l)*	0.135/0.133	
1G-H	Standard 3 (0.15 µg/l)*	0.319/0.329	
2A-B	Standard 4 (0.5 µg/l)*	0.907/0.880	
2C-D	Standard 5 (1.5 µg/l)*	2.405/2.304	
2E-F	Unknown 1	0.184/0.185	0.079
2G-H	Unknown 2	0.571/0.560	0.290
3A-B	Unknown 3	1.597/1.543	0.977

\*Exact concentration indicated on vial label.

### LIMITATIONS OF THE PROCEDURE

Grossly lipemic, icteric or hemolyzed samples do not interfere in the assay.

### EXPECTED VALUES

Good practice dictates that each laboratory establishes its own expected range of values.

### PERFORMANCE CHARACTERISTICS

#### 1.3 Detection limit

Detection limit is defined as the Capability of Detection according to ISO11843-Part 1. Capability of Detection should be seen as part of a method validation, rather than the lowest concentration that can be measured.

The detection limit is 0.01 (µg/l) determined with the methodology described in ISO11843-Part 4. Concentration of samples with absorbance below Standard 1 should not be calculated, instead expressed as less or equal to ( $\leq$ ) the concentration indicated on the vial for Standard 1.

#### 1.4 Recovery

Recovery upon addition is 100 - 111 % (mean 105 %).

Recovery upon dilution is 90 - 112 % (mean 97 %).

## DRG<sup>®</sup> Insulin (Porcine) ELISA (EIA-4747)

As of 10 Sept. 2007

For Veterinary Use Only

### 1.5 Hook effect

Samples with a concentration of up to 1000 µg/l can be measured without giving falsely low results.

### 1.6 Precision

Each sample was analyzed in 4-replicates on 20 different occasions.

Sample	Mean value (µg/l)	Coefficient of variation		
		within assay %	between assay %	total assay %
1	0.062	3.5	3.9	4.3
2	0.127	3.8	2.6	3.2
3	0.239	3.1	1.5	2.2

### 1.7 Specificity

The following cross reactions have been found:

Porcine C-peptide	< 0.001 %
Porcine Proinsulin	< 0.2 %
Human Insulin	28%
Human C-peptide	< 0.01 %
Human Proinsulin	< 0.1 %
NovoRapid	0.7 %
Levemir	< 0.0000002 %
Lantus	5.4 %
Humalog	< 0.0000006 %

## CALIBRATION

The Insulin (Porcine) ELISA is calibrated against an inhouse reference preparation of porcine insulin.

## WARRANTY

The performance data presented here was obtained using the procedure indicated. Any change or modification in the procedure not recommended by DRG may affect the results, in which event DRG disclaims all warranties expressed, implied or statutory, including the implied warranty of merchantability and fitness for use. DRG and its authorized distributors, in such event, shall not be liable for damages indirect of consequential.

## DRG<sup>®</sup> Insulin (Porcine) ELISA (EIA-4747)

As of 10 Sept. 2007

For Veterinary Use Only

### REFERENCES



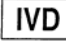
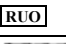

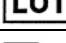
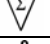






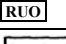

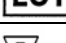
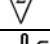



1. Bogdani M, Suenens K, Bock T, Pipeleers-Marichal M, In't Veld P, Pipeleers D (2005) Growth and Functional Maturation of beta-cells in Implants of Endocrine Cells Purified From Prenatal Porcine Pancreas. *Diabetes* 54:3387-3394
2. Dufrane D, D'hoore W, Goebbles RM, Saliez A, Guiot Y, Gianello P (2006) Parameters favouring successful adult pig islet isolations for xenotransplantation in pig-to-primate models. *Xenotransplantation* 13:204-214
3. Groth CG, Tibell A, Wennberg L, Bennet W, Lundgren T, Rydgård KJ, Lundin S, Lindeborg E, Korsgren O (2000) Clinical aspects and perspectives in islet xenotransplantation. *J Hepatobiliary Pancreat Surg* 7:364-369
4. Korsgren O, Jansson L, Eizirik D, Andersson A (1991) Functional and morphological differentiation of fetal porcine islet-like cell clusters after transplantation into nude mice. *Diabetologia* 34:379-386
5. Martinez-Puig D, Mourot J, Ferchaud-Roucher V, Anguita M, Garcia F, Krempf M, Pérez, JF (2006) Consumption of resistant starch decreases lipogenesis in adipose tissues but not in muscular tissues of growing pigs. *Livestock Science* 99:237-247
6. Ramsay TG, Bush JA, McMurtry JP, Thivierge MC, Davis TA (2004) Peripheral leptin administration alters hormone and metabolite levels in the young pig. *Comp Biochem Physiol* 138:17-25

# DRG® Insulin (Porcine) ELISA (EIA-4747)

As of 10 Sept. 2007

For Veterinary Use Only

## SYMBOLS USED WITH DRG ASSAY'S

Symbol	English	Deutsch	Francais	Español	Italiano
	Consult instructions for use	Gebrauchsanweisung beachten	Consulter les instructions d'utilisation	Consulte las instrucciones de uso	Consultare le istruzioni per l'uso
	European Conformity	CE-Konformitätskennzeichnung	Conformité aux normes européennes	Conformidad europea	Conformità europea
	In vitro diagnostic device	In-vitro-Diagnostikum	Usage Diagnostique in vitro	Para uso Diagnóstico in vitro	Per uso Diagnostica in vitro
	For research use only	Nur für Forschungszwecke	Seulement dans le cadre de recherches	Sólo para uso en investigación	Solo a scopo di ricerca
	Catalogue number	Katalog-Nr.	Numéro de catalogue	Número de catálogo	Numero di Catalogo
	Lot. No. / Batch code	Chargen-Nr.	Numéro de lot	Número de lote	Numero di lotto
	Contains sufficient for <n> tests/	Ausreichend für "n" Ansätze	Contenu suffisant pour "n" tests	Contenido suficiente para <n> ensayos	Contenuto sufficiente per "n" saggi
	Storage Temperature	Lagerungstemperatur	Temperature de conservation	Temperatura de conservación	Temperatura di conservazione
	Expiration Date	Mindesthaltbarkeits-datum	Date limite d'utilisation	Fecha de caducidad	Data di scadenza
	Legal Manufacturer	Hersteller	Fabricant	Fabricante	Fabbricante
Distributed by	Distributor	Distributeur	Distributeur	Distribuidor	Distributore
Content	Content	Inhalt	Conditionnement	Contenido	Contenuto
Volume/No.	Volume / No.	Volumen/Anzahl	Volume/Quantité	Volumen/Número	Volume/Quantità
Symbol	Portugues	Dansk	Svenska	Ελληνικά	
	Consulte as instruções de utilização	Se brugsanvisning	Se bruksanvisningen	Εγχειρίδιο χρήστη	
	Conformidade com as normas europeias	Europaeisk overensstemmelse	Europeisk överensstämmelse	Ευρωπαϊκή Συμμόρφωση	
	Diagnóstico in vitro	In vitro diagnostik	Diagnostik in vitro	in vitro διαγνωστικό	
					
	Catálogo n.º	Katalognummer	Katalog nummer	Αριθμός καταλόγου	
	No do lote	Lot nummer	Batch-nummer	Αριθμός Παρτίδος	
		Indeholder tilstrækkeligt til "n" test	Innehåller tillräckligt till "n" tester	Περιεχόμενο επαρκές για «n» εξετάσεις	
	Temperatura de conservação	Opbevarings-temperatur	Förvaringstemperatur	Θερμοκρασία αποθήκευσης	
	Prazo de validade	Udløbsdato	Bäst före datum	Ημερομηνία λήξης	
	Fabricante	Producent	Tillverkare	Κατασκευαστής	
Distributed by					
Content	Conteúdo	Indhold	Innehåll	Περιεχόμενο	
Volume/No.	Volume/Número	Volumen/antal	Volym/antal	Όγκος/αριθ.	