

DRG® Met Combi (Urine) ELISA (EIA-4083)

Revised 20 Apr. 2010 rm (Vers. 8.0)

RUO in the USA

Please use only the valid version of the package insert provided with the kit.

1 INTENDED USE AND PRINCIPLE OF THE TEST

Enzyme Immunoassay for the quantitative determination of Metanephine and Normetanephine in urine.

First Metanephine (Metadrenaline) and Normetanephine (Normetadrenaline) are quantitatively acylated.

The subsequent competitive ELISA kit uses the microtiter plate format. The antigen is bound to the solid phase of the microtiter plate. The acylated standards, controls and samples and the solid phase bound analytes compete for a fixed number of antiserum binding sites. After the system is in equilibrium, free antigen and free antigen-antiserum complexes are removed by washing. The antibody bound to the solid phase is detected by an anti-rabbit IgG-peroxidase conjugate using TMB as a substrate. The reaction is monitored at 450 nm.

Quantification of unknown samples is achieved by comparing their absorbance with a reference curve prepared with known standard concentrations.

During the sample preparation Metanephine (Metadrenaline) is quantitatively acylated.

2 ADVICE ON HANDLING THE TEST

2.1 Reliability of the test results

In order to assure a reliable evaluation of the test results it must be conducted according to the instructions included and in accordance with current rules and guidelines (GLP, RILIBÄK, etc.). Special attention must be paid to control checks for precision and correctness during the test; the results of these control checks have to be within the norm range. In case of significant discrepancies between the pre-set assay characteristics of this test and the actual results please contact the manufacturer of the test kit for further instructions.

It is recommended that each laboratory establishes its own reference intervals. The values reported in this test instruction are only indicative.

The results obtained with this test kit should not be taken as the sole reason for any therapeutic consequence but have to be correlated to other diagnostic tests and clinical observations.

2.2 Complaints

In case of complaints please submit to the manufacturer a written report containing all data as to how the test was conducted, the results received and a copy of the original test printout. Please contact the manufacturer to obtain a reclamation form and return it completely filled in to the manufacturer.

2.3 Warranty

This test kit was produced according to the latest developments in technology and subjected to stringent internal and external quality control checks. Any alteration of the test kit or the test procedure as well as the usage of reagents from different charges may have a negative influence on the test results and are therefore not covered by warranty. The manufacturer is not liable for damages incurred in transit.

2.4 Disposal

Residual substances and/or all remaining chemicals, reagents and ready for use solutions, are special refuse. The disposal is subject to the laws and regulations of the federation and the countries. About the removal of special refuse the responsible authorities or refuse disposal enterprises inform. The disposal of the kit must be made according to the national official regulations. Legal basis for the disposal of special refuse is the cycle economic- and waste law.

DRG[®] Met Combi (Urine) ELISA (EIA-4083)

Revised 20 Apr. 2010 rm (Vers. 8.0)

RUO in the USA

The appropriate safety data sheets of the individual products are available upon request. The safety data sheets correspond to the standard: ISO 11014-1.

2.5 Interference

Do not mix reagents and solutions from different lots. Consider different transport and storage conditions. Inappropriate handling of test samples or deviations from the test regulation can the results affect. Use no kit components beyond the expiration date. Avoid microbiological contamination of the reagents and the washing water. Consider incubation periods and wash references.

2.6 Precautions

Observe the incubation periods and washing instructions. Never pipette by mouth and avoid contact of reagents and specimens with skin. No smoking, eating or drinking in areas where samples or kit test tubes are handled. When working with kit components or samples, always wear protective gloves and wash your hand thoroughly as soon as you have finished the work. Avoid spraying of any kind. Avoid any skin contact with reagents. Use protective clothing and disposable gloves.

All steps have to be performed according to the protocol. Optimal test results are only obtained when using calibrated pipettes. Sodium azide could react with lead and copper tubes and may form highly explosive metal azide. When clearing up, rinse thoroughly with large volumes of water to prevent such formation.

All reagents of this testkit which contain human or animal serum or plasma have been tested and confirmed negative for HIV I/II, HbsAg and HCV by FDA approved procedures.

All reagents, however, should be treated as potential biohazards in use and for disposal.

3 STORAGE AND STABILITY

Store the reagents at 2 - 8 °C until expiration date. Do not use components beyond the expiry date indicated on the kit labels. Do not mix various lots of any kit component within an individual assay.

DRG® Met Combi (Urine) ELISA (EIA-4083)

Revised 20 Apr. 2010 rm (Vers. 8.0)

RUO in the USA

4 CONTENTS OF THE KIT

REAC-TUBES	Reaction Tubes	2 x 50	ready for use
WASH-CONC 50x	Wash Buffer Concentrate	1 x 20 mL	concentrate, dilute content with dist. water to a final volume of 1000 mL
CONJUGATE	Enzyme Conjugate	2 x 12 mL	ready for use, anti-rabbit IgG conjugated with peroxidase
SUBSTRATE	Substrate	2 x 12 mL	ready for use, containing a solution of tetramethylbenzidine (TMB)
STOP-SOLN	Stop Solution	2 x 12 mL	ready for use, containing 0.25 M H ₂ SO ₄
ARD MN	Adrenaline-Metanephrine Microtiter Strips	1 x 96 wells	12 strips, 8 wells each, break apart, pre-coated, blue coloured
NAD NMN	Noradrenaline-Normetanephrine Microtiter Strips	1 x 96 wells	12 strips, 8 wells each, break apart, pre-coated, yellow coloured
NMN-AS	Normetanephrine Antiserum	1 x 12 mL	from rabbit, ready for use, yellow coloured, yellow screw cap
MN-AS	Metanephrine Antiserum	1 x 12 mL	from rabbit, ready for use, blue coloured, blue screw cap
ACYL-CONC	Acylation Concentrate	1 x 0.5 mL	Concentrate. Has to be diluted prior to use.
ACYL-DILUENT	Acylation Diluent	1 x 4 mL	ready for use
STANDARD A	Standard A	1 x 2 mL	ready for use
STANDARD B	Standard B	1 x 2 mL	ready for use
STANDARD C	Standard C	1 x 2 mL	ready for use
STANDARD D	Standard D	1 x 2 mL	ready for use
STANDARD E	Standard E	1 x 2 mL	ready for use
STANDARD F	Standard F	1 x 2 mL	ready for use
HCL	Hydrochloric Acid	1 x 30 mL	ready for use, contains 0.25 M HCl
CONTROL 1	Control 1	1 x 2 mL	ready for use
CONTROL 2	Control 2	1 x 2 mL	ready for use
ACYL-BUFF	Acylation Buffer	1 x 30 mL	ready for use

4.1 Additional materials and equipment required but not provided with the kit

- Calibrated variable precision micropipettes (e.g. 10-100 µL / 100-1.000µL)
- Microtiter plate washing device
- ELISA reader capable of reading absorbance at 450 nm and 620 or 650 nm
- Centrifuge capable of at least 3.000 x g
- Absorbent material (paper towel)
- Distilled water, Vortex mixer, Temperature controlled water bath (90°C) or similar heating device

DRG International, Inc., USA Fax: (908) 233 0758 e-mail: corp@drg-international.com

Website: www.drg-international.com

DRG[®] Met Combi (Urine) ELISA (EIA-4083)

Revised 20 Apr. 2010 rm (Vers. 8.0)

RUO in the USA

■ The assay can be performed with or without shaking. If a shaker is used, it should have the following characteristics: shaking amplitude 3mm; approx. 600 rpm

5 SAMPLE COLLECTION AND STORAGE

Spontaneous or 24-hour urine, collected in a bottle containing 10-15 mL of 6 M HCl, should be used.

Determine the total volume of urine excreted during a period of 24 h for calculation of the results.

Storage: for longer periods (up to 6 months) at -20°C.

Repeated freezing and thawing should be avoided.

Avoid exposure to direct sunlight.

6 TEST PROCEDURE

Allow all reagents to reach room temperature and mix thoroughly by gentle inversion before use. Number the Reaction Tubes accordingly. Duplicate determinations are recommended.

■ *The sample preparation (hydrolysis and acylation) is identical for both the Metanephrine and Normetanephrine assay and has to be done only once.*

6.1 Preparation of reagents**Wash Buffer**

Dilute the 20 mL Wash Buffer Concentrate with distilled water to a final volume of 1000 mL.

Storage: up to 6 months 2–8°C

Acylation Solution

■ *Before preparing the Acylation Solution make sure that the Acylation Diluent has reached room temperature ($\geq 20^{\circ}\text{C}$) and forms a homogenous, crystal-free solution.*

Dilute the Acylation Concentrate 1 + 60 with Acylation-Diluent in a glass or polypropylene-vial.

Acylation Concentrate	10 μL	20 μL	25 μL	50 μL
Acylation-Diluent	600 μL	1.2 mL	1.5 mL	3 mL

■ *The Acylation Solution has to be prepared freshly prior to the assay (not longer than 60 minutes in advance). Discard after use!*

6.2 Sample preparation and acylation**Hydrolysis**

1. Pipette **25 μL** of standards, **25 μL** of controls, and **25 μL** of urine samples into the respective **Reaction Tubes**.
2. Add **250 μL** Hydrochloric Acid to all tubes.
3. Mix thoroughly (vortex) and hydrolyze for **30 min.** at **90 °C**.

DRG® Met Combi (Urine) ELISA (EIA-4083)

Revised 20 Apr. 2010 rm (Vers. 8.0)

RUO in the USA

4. Cool down the tubes to room temperature.

▲ *For the measurement of the free Metanephrine and free Normetanephrine only, leave away step 3 and 4.*

Acylation

1. Pipette **250 µL** of **Acylation Buffer** into all tubes.
2. Add **25 µL** of **Acylation Solution** to all tubes.
3. Mix thoroughly (vortex) and acylate for **15 minutes** at **RT** (20-25°C).
4. Add **2.5 mL dist. water** to all tubes.

▲ Take **25 µL** of the **acylated standards, controls and urine samples** for the **Metanephrine ELISA** and **Normetanephrine ELISA**.

6.3 Metanephrine ELISA

The usage of a shaker is not mandatory. The alternative protocol without shaker is highlighted in italic and shaded in grey.

1. Pipette **25 µL** of the **acylated standards, controls and samples** into the appropriate wells of the **Metanephrine Microtiter Strips**.
2. Pipette **100 µL** of the **Metanephrine Antiserum** into all wells.
3. Incubate **30 min** at **RT** (20-25°C) on a shaker (approx. 600 rpm).

Without usage of a shaker: shake the Metanephrine Microtiter Strips shortly by hand and incubate for 1 hour at RT (20-25°C).

4. Discard or aspirate the contents of the wells and **wash** each well **3 times** thoroughly with **300 µL Washbuffer**. Blot dry by tapping the inverted plate on absorbent material.
5. Pipette **100 µL** of the **Enzyme Conjugate** into all wells.
6. Incubate for **15 min** at **RT** (20-25°C) on a shaker (approx. 600 rpm).

Without usage of a shaker: incubate for 15 min at RT (20-25°C).

7. Discard or aspirate the contents of the wells and **wash** each well **3 times** thoroughly with **300 µL Washbuffer**. Blot dry by tapping the inverted plate on absorbent material.
8. Pipette **100 µL** of the **Substrate** into all wells.
9. Incubate for **15 ± 2 min** at **RT** (20-25°C) on a shaker (approx. 600 rpm).

▲ ***Without usage of a shaker: incubate for 15 min ± 2 at RT (20-25°C).***

Avoid exposure to direct sun light!

10. Add **100 µL** of the **Stop Solution** to each well and shake the microtiter plate to ensure a homogeneous distribution of the solution.
11. **Read** the absorbance of the solution in the wells within 10 minutes, using a microplate reader set to **450 nm** and a reference wavelength between 620 nm and 650 nm.

6.4 Normetanephrine ELISA

The usage of a shaker is not mandatory. The alternative protocol without shaker is highlighted in italic and shaded in grey.

DRG® Met Combi (Urine) ELISA (EIA-4083)

Revised 20 Apr. 2010 rm (Vers. 8.0)

RUO in the USA

1. Pipette **25 µL** of the **acylated standards, controls and samples** into the appropriate wells of the **Normetanephine Microtiter Strips**.
2. Pipette **100 µL** of the **Normetanephine Antiserum** into all wells.
3. Incubate **30 min at RT** (20-25°C) on a shaker (approx. 600 rpm).
Without usage of a shaker: shake the Normetanephine Microtiter Strips shortly by hand and incubate for 1 hour at RT (20-25°C).
4. Discard or aspirate the contents of the wells and **wash** each well **3 times** thoroughly with **300 µL Washbuffer**. Blot dry by tapping the inverted plate on absorbent material.
5. Pipette **100 µL** of the **Enzyme Conjugate** into all wells.
6. Incubate for **15 min at RT** (20-25°C) on a shaker (approx. 600 rpm).
Without usage of a shaker: incubate for 15 min at RT (20-25°C).
7. Discard or aspirate the contents of the wells and **wash** each well **3 times** thoroughly with **300 µL Washbuffer**. Blot dry by tapping the inverted plate on absorbent material.
8. Pipette **100 µL** of the **Substrate** into all wells.
9. Incubate for **15 ± 2 min at RT** (20-25°C) on a shaker (approx. 600 rpm).
Without usage of a shaker: incubate for 15 min ± 2 at RT (20-25°C).
Avoid exposure to direct sun light!
10. Add **100 µL** of the **Stop Solution** to each well and shake the microtiter plate to ensure a homogeneous distribution of the solution.
11. **Read** the absorbance of the solution in the wells within 10 minutes, using a microplate reader set to **450 nm** and a reference wavelength between 620 nm and 650 nm.

7 CALCULATION OF RESULTS

Standard	Concentration of the standards					
	A	B	C	D	E	F
Normetanephine (ng/mL=µg/L)	0	30	90	300	900	3 000
Normetanephine (nmol/L)	0	164	491	1 638	4 914	16 380
Metanephine (ng/mL=µg/L)	0	20	60	200	600	2 000
Metanephine (nmol/L)	0	101	304	1 014	3 042	10 140
Conversion:	Normetanephine (ng/mL) x 5.46 = Normetanephine (nmol/L) Metanephine (ng/mL) x 5.07 = Metanephine (nmol/L)					

The calibration curve is obtained by plotting the absorbance readings (calculate the mean absorbance) of the standards (linear, y-axis) against the corresponding standard concentrations (logarithmic, x-axis).

Use a non-linear regression for curve fitting (e.g. spline, 4- parameter, akima).

DRG[®] Met Combi (Urine) ELISA (EIA-4083)

Revised 20 Apr. 2010 rm (Vers. 8.0)

RUO in the USA

The concentrations of the samples can be read directly from the standard curve.

The amount of analyte excreted per day (µg/day) is calculated according to:

concentration of the sample (in µg/L) x volume of urine excreted per day (in L/day)

Example

The concentration of the sample read from the curve is 125 µg/L. The amount of urine collected during 24 hours is 1.3 L. Then the amount of analyte excreted during one day would be:


$$125 \mu\text{g/L} \times 1.3 \text{ L/day} = 162.5 \mu\text{g/day}$$

7.1 Quality control

It is recommended to use control samples according to state and federal regulations. Use controls at both normal and pathological levels. The kit or other commercial controls should fall within established confidence limits. The confidence limits of the kit controls are indicated on the QC Report.

7.2 Calibration

The binding of the antisera and the enzyme conjugates and the activity of the enzyme used are temperature dependent, and the extinction values may vary if a thermostat is not used. The higher the temperature, the higher the extinction values will be. The extinction values also depend on the incubation times. The optimal temperature during the Enzyme Immunoassay is between 20-25°C.



 *In case of overflow, read the absorbance of the solution in the wells within 10 minutes, using a microplate reader set to 405 nm*

DRG[®] Met Combi (Urine) ELISA (EIA-4083)

Revised 20 Apr. 2010 rm (Vers. 8.0)

RUO in the USA

7.3 Typical calibration curves

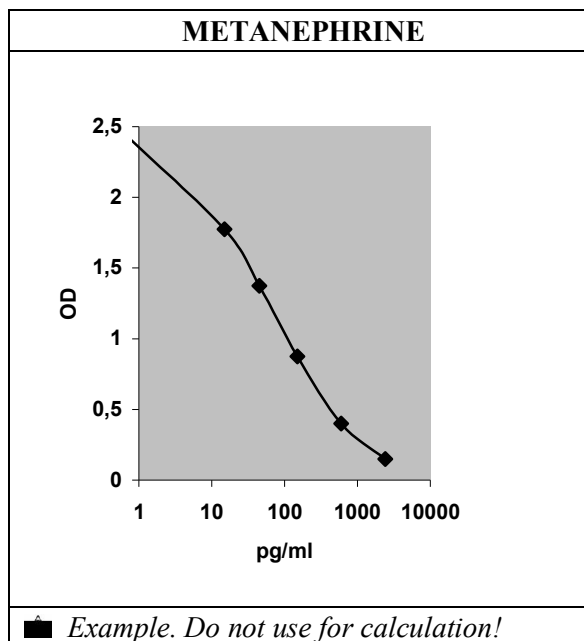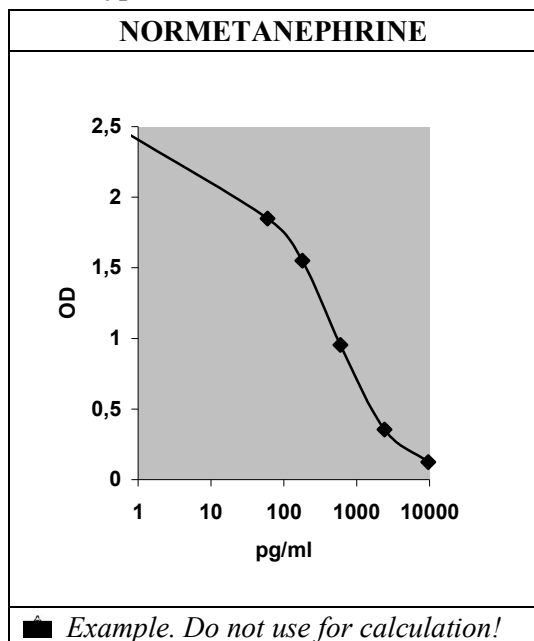

8 ASSAY CHARACTERISTICS

Expected Reference Values		Metanephine	Normetanephine
	Urine	< 350 µg/day	< 600 µg/day

Analytical Sensitivity (Limit of Detection)		Metanephine	Normetanephine
	Urine	13 ng/mL	23 ng/mL

Analytical Specificity (Cross Reactivity)	Substance	Cross Reactivity (%)	
		Metanephine	Normetanephine
	Derivatized Metanephine	100	0.11
	Derivatized Normetanephine	0.15	100
	Derivatized 3-methoxytyramine	< 0.01	0.19
	Adrenaline	3.3	< 0.001
	Noradrenaline	< 0.001	0.64
	Dopamine	< 0.001	< 0.01
	Vanillic mandelic acid, L-Dopa, Homovanillic acid, L-Tyrosin, Tyramin	< 0.001	< 0.001

DRG® Met Combi (Urine) ELISA (EIA-4083)

Revised 20 Apr. 2010 rm (Vers. 8.0)

RUO in the USA


Precision							
Intra-Assay				Inter-Assay			
	Sample	Range (ng/mL)	CV (%)		Sample	Range (ng/mL)	CV (%)
Metanephrine	1	69 ± 8.6	12.6	Metanephrine	1	102 ± 15.4	15.1
	2	446 ± 23	5.2		2	448 ± 40	8.9
Normetanephrine	1	200 ± 34	17.2	Normetanephrine	1	191 ± 41	21.4
	2	857 ± 153	17.8		2	682 ± 131	19.3

Linearity			Range	Serial dilution up to	Mean (%)
	Metanephrine	Urine	40 – 1.600 ng/mL	1:16	98
	Normetanephrine	Urine	40 – 5 200 ng/mL	1:16	93

Recovery			Mean (%)	Range (%)	% Recovery after spiking
	Metanephrine	Urine	105	98 – 119	
	Normetanephrine	Urine	103	90 – 113	

Method Comparison versus HPLC*	Metanephrine	Urine	HPLC = 0.9 ELISA – 0.8	r = 0.99; n = 40
	Normetanephrine	Urine	HPLC = 0.9 ELISA + 0.6	r = 0.99; n = 40

* The concentrations were assessed using both the ELISA and the HPLC method (external QC samples from UK NEQAS). The correlation between ELISA and HPLC is excellent. Please take in mind, that the UK control values are the mean of about 40 different HPLC users, and contain always one pathological sample per sending.




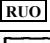






 **For actual literature, information about clinical significance or any other information please contact your local supplier.**





DRG® Met Combi (Urine) ELISA (EIA-4083)

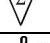Revised 20 Apr. 2010 rm (Vers. 8.0)

RUO in the USA

SYMBOLS USED WITH DRG ASSAYS

Symbol	English	Deutsch	Français	Español	Italiano
	Consult instructions for use	Gebrauchsanweisung beachten	Consulter les instructions d'utilisation	Consulte las instrucciones de uso	Consultare le istruzioni per l'uso
	European Conformity	CE-Konformitätskennzeichnung	Conformité aux normes européennes	Conformidad europea	Conformità europea
	In vitro diagnostic device	In-vitro-Diagnostikum	Usage Diagnostic in vitro	Para uso Diagnóstico in vitro	Per uso Diagnostica in vitro
	For research use only	Nur für Forschungszwecke	Seulement dans le cadre de recherches	Sólo para uso en investigación	Solo a scopo di ricerca
	Catalogue number	Katalog-Nr.	Numéro de catalogue	Número de catálogo	Numero di Catalogo
	Lot. No. / Batch code	Chargen-Nr.	Numéro de lot	Número de lote	Numero di lotto
	Contains sufficient for <n> tests/	Ausreichend für "n" Ansätze	Contenu suffisant pour "n" tests	Contenido suficiente para <n> ensayos	Contenuto sufficiente per "n" saggi
	Storage Temperature	Lagerungstemperatur	Température de conservation	Temperatura de conservación	Temperatura di conservazione
	Expiration Date	Mindesthaltbarkeits-datum	Date limite d'utilisation	Fecha de caducidad	Data di scadenza
	Legal Manufacturer	Hersteller	Fabricant	Fabricante	Fabbricante
Distributed by	Distributor	Vertreiber	Distributeur	Distribuidor	Distributore
Content	Content	Inhalt	Conditionnement	Contenido	Contenuto
Volume/No.	Volume / No.	Volumen/Anzahl	Volume/Quantité	Volumen/Número	Volume/Quantità

Symbol	Portugues	Dansk	Svenska	Ελληνικά
	Consulte as instruções de utilização	Se brugsanvisning	Se bruksanvisningen	Εγχειρίδιο χρήστη
	Conformidade com as normas europeias	Europæisk overensstemmelse	Europeisk överensstämmelse	Ευρωπαϊκή Συμμόρφωση
	Diagnóstico in vitro	In vitro diagnostik	Diagnostik in vitro	in vitro διαγνωστικό
	Catálogo n.º	Katalognummer	Katalog nummer	Αριθμός καταλόγου
	No do lote	Lot nummer	Batch-nummer	Αριθμός Παρτίδος
		Indeholder tilstrækkeligt til "n" test	Innehåller tillräckligt till "n" tester	Περιεχόμενο επαρκές για «n» εξετάσεις
	Temperatura de conservação	Opbevarings-temperatur	Förvaringstemperatur	Θερμοκρασία αποθήκευσης
	Prazo de validade	Udløbsdato	Bäst före datum	Ημερομηνία λήξης
	Fabricante	Producent	Tillverkare	Κατασκευαστής
Content	Conteúdo	Indhold	Innehåll	Περιεχόμενο
Volume/No.	Volume/Número	Volumen/antal	Volym/antal	Όγκος/αριθ.