



Peninsula Laboratories, LLC

A Member of the Bachem Group

305 Old County Road, San Carlos, CA 94070

Tel: (800) 922-1516 • (650) 592-5392

Fax: (650) 595-4071

www.bachem.com



Cat. No. S-1282 EIA for D-Marker peptide (for human serum and plasma samples)

SPECIFICATIONS: Extraction-free kit for human serum/plasma


STANDARD: 1 µg (not available separately) H-D-Cys-Gly-D-Pro-Gly-Gly-D-Tyr-D-Lys-D-Val-D-Leu-D-Val-D-Leu-D-Asn-D-Pro-D-Ser-D-Val-D-Ala-D-Ala-D-Thr-OH

PROTOCOL: V (Ab1hr.Std2hr.BtON)

RANGE: 0-25 ng/ml

AVERAGE IC₅₀: 0.4 ng/ml

SHELF-LIFE: One year from production date

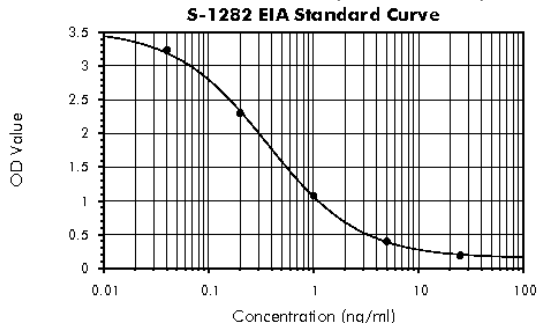

Suggested Preparation of Standards		
	ng/ml	Range: 0.02 to 25 ng/ml
Stock	1000	
S1	25.00	Add 25 µl stock + 975 µl diluent
S2	6.25	Add 200 µl S1 + 600 µl diluent
S3	1.56	Add 200 µl S2 + 600 µl diluent
S4	0.39	Add 200 µl S3 + 600 µl diluent
S5	0.10	Add 200 µl S4 + 600 µl diluent
S6	0.02	Add 200 µl S5 + 600 µl diluent
S0	0.00	500 µl diluent

CROSS REACTION:

PEPTIDE

D-Marker peptide

%

100