

Code No. 27415

Mouse/Rat sAPP α Assay Kit - IBL**INTRODUCTION**

Alzheimer's disease (AD) was first reported by A. Alzheimer, a German neuropathologist in 1907 and is considered as a major factor of dementia. It is known that Amyloid β (A β); which is major constituent of senile plaque) is cleaved from Amyloid Precursor Protein (APP; which exists in three main isoforms, APP695, APP751, and APP770) by β -secretase and subsequent γ -secretase (ref. 1). The production of soluble APP β (sAPP β) by β -secretase cleavage corresponds to A β production accordingly, so it is desired to measure sAPP β in parallel with A β . On the other hand, it is considered that in the metabolic pathway of APP, APP is first cleaved by α -secretase rather than β -secretase normally to produce soluble APP α (sAPP α) and subsequently P3 is cleaved from the remaining C-terminal fragment by γ -secretase. In recent research, there are several attempts to apply the inhibitor of β -secretase and the activation of α -secretase for AD treatment. As stated above, it is useful for research of AD to measure sAPP α as well as sAPP β .

This kit can measure Mouse/Rat sAPP α .

Note: Please pay attention in sample selection since this kit also measures full-length APP.

PRINCIPLE

This kit is a solid phase sandwich ELISA using 2 kinds of high specific antibodies. Tetra Methyl Benzidine (TMB) is used as a coloring agent (Chromogen). The strength of coloring is proportional to the quantities of Mouse/Rat sAPP α .

MEASUREMENT RANGE

10 - 640 pg/mL

INTENDED USE

- This IBL's assay kit is capable for the quantitative determination each of mouse and rat sAPP α in serum, EDTA plasma, cell culture supernatant and brain extract.
- Note: Please pay attention in measuring of brain samples since this kit also measures full length APP molecules.
- When the culture media samples contain serum like FCS, any cross-reaction may be observed. So we recommend you to set the negative control.

KIT COMPONENT

- | | | | |
|---|-------------------------------|--|------------|
| 1 | Precoated plate | : Anti- Mouse/Rat APP (597) Rabbit IgG Affinity Purify | 96Well x 1 |
| 2 | Labeled antibody Conc. | | |
| | | : (30X) HRP conjugated Anti-Mouse/Rat APP (18) Rabbit IgG Fab' Affinity Purify | 0.4mL x 1 |
| 3 | Standard | : Recombinant Mouse/Rat sAPP α | 0.5mL x 2 |
| 4 | EIA buffer | : 1% BSA, 0.05% Tween20 in PBS | 30mL x 1 |
| 5 | Solution for Labeled antibody | : 1% BSA, 0.05% Tween20 in PBS | 12mL x 1 |
| 6 | Chromogen | : TMB solution | 15mL x 1 |
| 7 | Stop solution | : 1N H ₂ SO ₄ | 12mL x 1 |
| 8 | Wash buffer Conc. | : (40X) 0.05% Tween20 in phosphate buffer | 50mL x 1 |

OPERATION MANUAL**1. Materials needed but not supplied**

- Plate reader (450nm)
- Graduated cylinder and beaker
- Refrigerator (as 4°C)
- Paper towel
- Washing bottle for precoated plate
- Disposable test tube for "2, Labeled antibody Conc." and "6, Chromogen"
- Micropipette and tip
- Deionized water
- Graph paper (log/log)
- Tube for dilution of Standard

2. Preparation


- 1) Preparation of wash buffer
"8, Wash buffer Conc." is a concentrated (40X) buffer. Adjust the temperature of "8, Washing buffer Conc." to room temperature and then, mix it gently and completely before use. Dilute 50 mL of "8, Wash buffer Conc." with 1,950 mL of deionized water and mix it. This is the wash buffer for use. This prepared wash buffer shall be stored in refrigerator and used within 2 weeks after dilution.
- 2) Preparation of Labeled antibody
"2, Labeled antibody Conc." is a concentrated (30X). Dilute "2, Labeled antibody Conc." with "5, Solution for Labeled antibody" in 30 times according to required quantity into a disposable test tube. Use this resulting solution as Labeled antibody.
Example)
In case you use one strip (8 well), the required quantity of Labeled antibody is 800 μ L. (Dilute 30 μ L of "2, Labeled antibody Conc." with 870 μ L of "5, Solution for Labeled antibody" and mix it. And use the resulting solution by 100 μ L in each well.)
This operation should be done just before the application of Labeled antibody.
The remaining "2, Labeled antibody Conc." should be stored at 4°C in firmly sealed vial.
- 3) Preparation of Standard
Put just 0.5 mL of deionized water into the vial of "3, Standard" and mix it gently and completely. This solution is 1,280 pg/mL Mouse/Rat sAPP α standard.
- 4) Dilution of Standard
Prepare 8 tubes for dilution of "3, Standard". Put 230 μ L each of "4, EIA buffer" into the tube.
Specify the following concentration of each tube.

Tube-1	640 pg/mL
Tube-2	320 pg/mL
Tube-3	160 pg/mL
Tube-4	80 pg/mL
Tube-5	40 pg/mL

Tube-6	20pg/mL
Tube-7	10 pg/mL
Tube-8	0 pg/mL (Test Sample Blank)

Put 230 μ L of Standard solution into tube-1 and mix it gently. Then, put 230 μ L of tube-1 mixture into tube-2. Dilute two times standard solution in series to set up 7 points of diluted standard between 640 pg/mL and 10 pg/mL. Tube-8 is the test sample blank as 0 pg/mL.

See following picture.

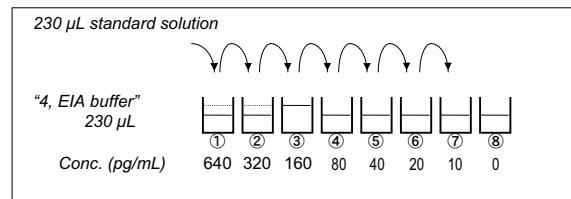**5) Dilution of test sample**

Test sample may be diluted with "4, EIA buffer" as necessary.

If the concentration of Mouse/Rat sAPP α in samples may not be estimated in advance, the pre-assay with several different dilutions will be recommended to determine the proper dilution of samples.

3. Measurement procedure

All reagents shall be brought to room temperature approximately 30 minutes before use. Then mix it gently and completely before use. Make sure of no change in quality of the reagents. Standard curve shall be prepared simultaneously with the measurement of test samples.

Reagents	Test Sample	Standard	Test Sample Blank	Reagent Blank
	Test sample 100 μ L	Diluted standard (Tube 1~7) 100 μ L	EIA buffer (Tube-8) 100 μ L	EIA buffer 100 μ L
Incubation overnight at 4°C with plate lid				
Washing 7 times				
Labeled Antibody	100 μ L	100 μ L	100 μ L	-
Incubation for 30 minutes at 4°C with plate lid				
Washing 9 times				
Chromogen	100 μ L	100 μ L	100 μ L	100 μ L
Incubation for 30 minutes at room temperature (shielded)				
Stop solution	100 μ L	100 μ L	100 μ L	100 μ L
Read the plate at 450nm against a Reagent Blank within 30 minutes after addition of Stop solution.				

- 1) Determine wells for reagent blank. Put 100 μ L each of "4, EIA buffer" into the wells.
- 2) Determine wells for test sample blank, test sample and diluted standard. Then, put 100 μ L each of test sample blank (tube-8), test sample and dilutions of standard (tube-1-7) into the appropriate wells.
- 3) Incubate the precoated plate overnight at 4°C after covering it with plate lid.
- 4) Wash each well of the precoated plate vigorously with wash buffer using the washing bottle. Then, fill each well with wash buffer and leave the precoated plate lay for 15-30 seconds. Remove wash buffer completely from the precoated plate by snapping. This procedure must be repeated more than 7 times. Then, remove the remaining liquid from all wells completely by snapping the precoated plate onto paper towel.
In case of using a plate washer, after 4 times washing with plate washer, washing with above washing bottle must be repeated 3 times.
- 5) Pipette 100 μ L of labeled antibody solution into the wells of test samples, diluted standard and test sample blank.
- 6) Incubate the precoated plate for 30 minutes at 4°C after covering it with plate lid.
- 7) Wash the precoated plate 9 times in the same manner as 4).
- 8) Taken the required quantity of "6, Chromogen" into a disposable test tube. Then, pipette 100 μ L from the test tube into the wells. Please do not return the rest of the test tube to "6, Chromogen" bottle to avoid contamination.
- 9) Incubate the precoated plate for 30 minutes at room temperature in the dark. The liquid will turn blue by addition of "6, Chromogen".
- 10) Pipette 100 μ L of "7, Stop solution" into the wells. Mix the liquid by tapping the side of precoated plate. The liquid will turn yellow by addition of "7, Stop solution".
- 11) Remove any dirt or drop of water on the bottom of the precoated plate and confirm there is no bubble on the surface of the liquid. Then, run the plate reader and conduct measurement at 450 nm against a reagent blank. The measurement shall be done within 30 minutes after addition of "7, Stop solution".

SPECIAL ATTENTION

- 1) Test samples should be measured soon after collection. For the storage of test samples, store them frozen and do not repeat freeze/thaw cycles. Thaw the test samples at a low temperature and mix them completely before measurement.
- 2) Test samples should be diluted with "4, EIA buffer", if the need arises.
- 3) Duplicate measurement of test samples and standard is recommended.
- 4) Use test samples in neutral pH range. The contaminations of organic solvent may affect the measurement.
- 5) Use only wash buffer contained in this kit for washing the precoated plate. Insufficient washing may lead to the failure in measurement.
- 6) Remove the wash buffer completely by tapping the precoated plate on paper


- towel. Do not wipe wells with paper towel.
- 7) "6, Chromogen" should be stored in the dark due to its sensitivity against light. "6, Chromogen" should be avoided contact with metals.
 - 8) Measurement should be done within 30 minutes after addition of "7, Stop solution".

CALCULATION OF TEST RESULT

Subtract the absorbance of test sample blank from all data, including standards and unknown samples before plotting. Plot the subtracted absorbance of the standards against the standard concentration on log-log graph paper. Draw the best smooth curve through these points to construct the standard curve. Read the concentration for unknown samples from the standard curve.

Example of standard curve

Conc. (pg/mL)	Absorbance (450nm)
640	2.077
320	1.109
160	0.587
80	0.324
40	0.186
20	0.120
10	0.087
0 (Test Sample Blank)	0.047

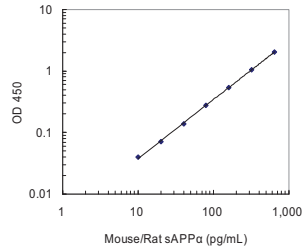

* The typical standard curve is shown above. This curve can not be used to derive test results. Please run a standard curve for each assay.

PERFORMANCE CHARACTERISTICS

1. Titer Assay (Samples with standard added are used.)

Specimen	Titer (x)	Measurement Value (pg/mL)	Theoretical Value (pg/mL)	%
10%FCS added RPMI-1640	4	145.93	160.00	91.2
	8	71.10	80.00	88.9
	16	36.65	40.00	91.6
Mouse Serum (BALB/c)	4	106.81	190.79	56.0
	8	51.35	90.73	56.6
	16	28.48	45.12	63.1
Mouse Plasma (EDTA) (BALB/c)	4	135.65	183.33	74.0
	8	66.56	86.02	77.4
	16	32.57	41.19	79.1
Mouse Brain Extract (BALB/c)	16	440.27	446.20	98.7
	32	161.84	168.20	96.2
	64	62.37	66.37	94.0
Rat Serum (SD)	2	71.16	338.22	21.0
	4	53.02	166.52	31.8
	8	35.35	82.93	42.6
Rat Plasma (EDTA) (SD)	2	297.38	349.34	85.1
	4	147.28	160.00	92.1
	8	73.60	80.00	92.0
Rat Brain Extract (SD)	16	257.22	275.09	93.5
	32	91.07	99.04	92.0
	64	40.76	44.96	90.7

2. Added Recovery Assay

Specimen	Theoretical Value (pg/mL)	Measurement Value (pg/mL)	%
10%FCS added RPMI-1640 (x4)	325.31	344.20	105.8
	165.31	175.22	106.0
	85.31	91.37	107.1
Mouse Serum (BALB/c) (x8)	169.29	104.36	61.6
	89.29	58.87	65.9
	49.29	34.57	70.1
Mouse Plasma (EDTA) (BALB/c) (x8)	171.01	142.07	83.1
	91.01	76.90	84.5
	51.01	45.38	89.0
Mouse Brain Extract (BALB/c) (x16)	431.28	465.03	107.8
	351.28	393.45	112.0
	311.28	328.99	105.7
Rat Serum (SD) (x4)	164.56	60.70	36.9
	84.56	35.13	41.5
	44.56	19.76	44.3

Specimen	Theoretical Value (ng/mL)	Measurement Value (ng/mL)	%
Rat Plasma (EDTA) (SD) (x4)	162.51	145.94	89.8
	82.51	74.08	89.8
	42.51	34.87	82.0
Rat Brain Extract (SD) (x16)	450.87	457.39	101.4
	370.87	369.82	99.7
	330.87	336.01	101.6

3. Intra - Assay

Measurement Value (pg/mL)	SD value	CV value (%)	n
402.49	11.73	2.9	24
112.47	4.76	4.2	24
31.81	2.84	8.9	24

4. Inter - Assay

Measurement Value (pg/mL)	SD value	CV value (%)	n
398.90	16.42	4.1	7
111.58	5.01	4.5	7
31.21	2.18	7.0	7

5. Specificity

Compound	Cross-reactivity
Mouse/Rat sAPPα	100 %
Human sAPPα	0.6 %
Human sAPPβ	≤0.1 %

6. Sensitivity

2.86 pg/mL

The sensitivity for this kit was determined using the guidelines under the National Committee for Clinical Laboratory Standards (NCCLS) Evaluation Protocols. (National Committee for Clinical Laboratory Standards Evaluation Protocols, SC1, (1989) Villanova, PA: NCCLS.)

PRECAUTION FOR INTENDED USE AND/OR HANDLING

1. All reagents should be stored at 2 - 8°C. All reagents shall be brought to room temperature approximately 30 minutes before use.
2. "3, Standard" is lyophilized products. Be careful to open this vial.
3. "7, Stop solution" is a strong acid substance. Therefore, be careful not to have your skin and clothes contact "7, Stop solution" and pay attention to the disposal of "7, Stop solution".
4. Dispose used materials after rinsing them with large quantity of water.
5. Precipitation may occur in "2, Labeled antibody Conc.", however, there is no problem in the performance.
6. Wash hands after handling reagents.
7. Do not mix the reagents with the reagents from a different lot or kit.
8. Do not use expired reagents.
9. This kit is for research purpose only. Do not use for clinical diagnosis.

STORAGE AND THE TERM OF VALIDITY

Storage Condition : 2 - 8 °C
The expiry date is specified on outer box.

REFERENCE

1. Selkoe DJ. Normal and abnormal biology of the beta-amyloid precursor protein. *Annu Rev Neurosci.* 1994;17:489-517.
2. Citron M, Oltersdorf T, Haass C, McConlogue L, Hung AY, Seubert P, Vigo-Pelfrey C, Lieberburg I, Selkoe DJ. Mutation of the beta-amyloid precursor protein in familial Alzheimer's disease increases beta-protein production. *Nature.* 1992 Dec 17;360(6405):672-4.
3. Hsiao K, Chapman P, Nilsen S, Eckman C, Harigaya Y, Younkin S, Yang F, Cole G. Correlative memory deficits, Aβeta elevation, and amyloid plaques in transgenic mice. *Science.* 1996 Oct 4;274(5284):99-102.

Version. 1.3

IBL Amyloidβ, APP related Product Lines:

Code No.	Name	Volume
27711	Human Amyloidβ (1-42) Assay Kit - IBL	96 Well
27712	Human Amyloidβ (1-42) (N) Assay Kit - IBL	96 Well
27713	Human Amyloidβ (1-40) Assay Kit - IBL	96 Well
27714	Human Amyloidβ (1-40) (N) Assay Kit - IBL	96 Well
27718	Human Amyloidβ (1-40) (FL) Assay Kit - IBL	96 Well
27720	Mouse/Rat Amyloidβ (1-40) High Specific Assay Kit - IBL	96 Well
27716	Human Amyloidβ (N3pE-42) Assay Kit - IBL	96 Well
27418	Human Amyloidβ (N3pE-40) Assay Kit - IBL	96 Well
27729	Human Amyloidβ (1-x) Assay Kit - IBL	96 Well
27724	Human sAPPα Assay Kit - IBL	96 Well
27415	Mouse/Rat sAPPα Assay Kit - IBL	96 Well
27722	Human sAPPβ-Wild Type Assay Kit - IBL	96 Well
27723	Human sAPPβ-Swedish Type Assay Kit - IBL	96 Well
27776	Human APP βCTF Assay Kit - IBL	96 Well