



cDNA Library, Planaria

02-709 500 ng

This cDNA library (plasmid DNA) is constructed from *Planaria*-derived poly(A)⁺ RNA by the Linker-Primer method (Ref.1) by Professor Hiroshi Nojima of Research Institute for Microbial Diseases, Osaka University. This library is unidirectionally cloned by using the oligo (dT)₁₈ linker primer which contains the restriction enzyme site of *Not* I, and *Bam*HI (*Bgl* II)-*Sma* I adaptor.

The pAP3*neo* vector used in this library can express in mammalian cells as it contains SV40 promoter. It also contains Ori for the plasmid replication and f1 ori which is necessary for ssDNA synthesis in *E. coli*, and bacteriophage T7 and T3 promoter for RNA synthesis (see Figure). GenBank Accession No. <u>AB003468</u>

Application

PCR screening of known or unknown gene: Prepare the primers for the known or unknown gene (cDNA) and amplify the gene by PCR from this library followed by cloning to an appropriate vector. It is useful for large-scale protein productions, and preparation of probes, etc.

Standard amplifying conditions: 35 cycles of PCR reactions using 10-100 ng of cDNA as a template. (Change the quantity of template and the number of cycles depending on the expression levelof mRNA of the particular gene.)

Specification

Quantity: 500 ng (40 ng/ul, 13ul) in 10 mM Tris-HCl-1mM EDTA (pH 7.5)

Quality: 1) Number of independent clones: 14 x 106

2) Average insert size: longer than 1 kb

Storage: -20

References

- 1. Kobori, M., Ikeda, Y., Nara, H., Kato, M., Kumegawa, M., Nojima, H., and Kawashima, H. "Large scale isolation of osteoclast-specific genes by an improved method involving the preparation of a subtracted cDNA library." *Genes Cells* **3**: 459-475 (1998) PMID: 9753427
- 2. Sambrook, J. and Russell, DW. *Molecular Cloning* Chapter 11 "Preparation of cDNA libraries and gene identification." CSHL Press (2001)


Note

- * This library is to be used only by the purchaser. It is not allowed to amplify and transfer the library to a third person.
- * Related products: human tissue specific cDNA libraries and cDNA libraries of model organisms (See <u>HP</u>).



SCETI

Multiple cloning site (MCS) of pAP3neo

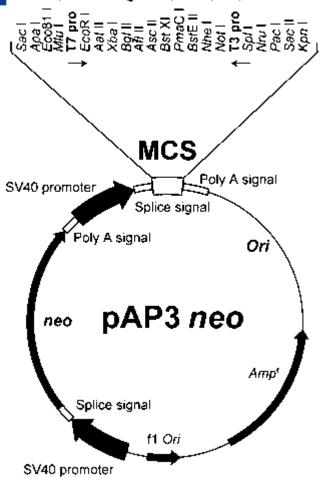


Multiple cloning site (MCS) of pAP3neo

Sacl Agal TTul TTul TTul Topo Agal Agal Notl Notl Notl Notl Notl Notl

Fig. Structure of pAP3neo and the restriction sites

Ori is an origin required for replication in *E. coli*

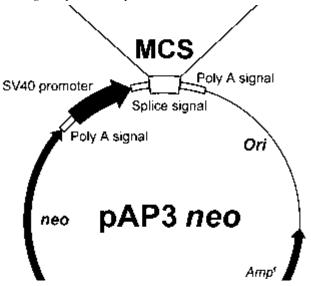



3-6-7 Kasumigaseki, Chiyoda-ku, Tokyo 100-0013, JAPAN Tel +81(3) 5510-2347 Fax +81(3) 5510-0133 URL: http://www.sceti.jp/export/ e-mail: exp-pet@sceti.co.jp