

Anti-Streptolysin O antibody, rabbit serum

Polyclonal antibody 64-001 50 μ l 64-002 250 μ l

Streptolysin O (SL0) is a membrane-damaging extracellular toxin produced by hemolytic streptococci. In the hemolytic streptococci infection, the antibodies to various proteins secreted from the bacteria such as toxins increase. For the diagnostic confirmation of SL0 production by the infection, the measurement of the antibody, ASLO (Anti-Streptolysin O), is most widely used.

The product was prepared by immunizing rabbit with purified recombinant SL0 of Group C hemolytic streptococci produced in *E. coli* as antigen. Using this antiserum it is possible to detect the SL0 secreted into the culture supernatant not only of Group A hemolytic streptolysin O but also of Group C and G origins.

Applications

- 1) Western blotting (x2,000~10,000 dilution)
- 2) Neutralization of hemolysis by Streptolysin O (SL0)
- 3) The anti-streptolysin O antibodies (ASLO) sold in the current market is obtained from SL0 antigen purified from the culture supernatant of hemolytic streptococci. But the antigens are contaminated with other secreted proteins. The present product was prepared from SL0 which was over-expressed in *E. coli* as recombinant protein and purified, thus showing extremely high specification (below figure).

Specifications

Reactivity: Streptolysin O (Streptolysin O of Group A, C, and G origins)

Form: Undiluted serum added with 0.09 % sodium azide added

Storage: 4°C

References:

- 1..H. Kimoto et al.: Expression of recombinant streptolysin O and specific antibody production. *J.Mol.Microbiol. Biotechnol.*,10,64-68 (2005).


Detection of streptolysin O secreted by hemolytic streptococci Group A in the culture supernatant by Western blotting , using anti-streptolysin O antibody.

Lane 1: Culture supernatant of hemolytic streptococci Group A

Lane 2: Culture media only (Control)

*This product is for research use only, not for use in human

Related product: 01-531 Streptolysin O

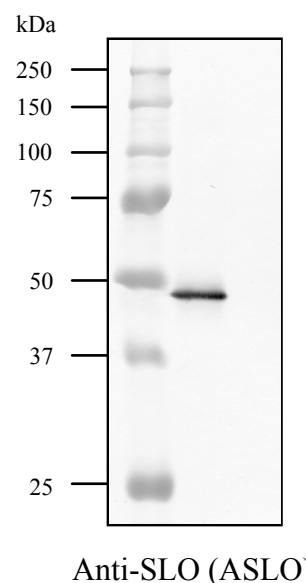

<Distributed by >: SCETI K.K. DF Kasumigaseki Place, 3-6-7 Kasumigaseki,
Chiyoda-ku Tokyo 100-0013 JAPAN Tel: +81-3-5510-2347 Fax: +81-3-5510-0134
E-mail: exp-pet@sceti.co.jp URL: www.sceti.co.jp/export/

<Manufactured by> : BioAcademia,Inc. 7-7-18 Saito-Asagi, Ibaraki, Osaka 567-0085, JAPAN